

EPX-3500 HD

Electronic Video Endoscopy System

VP-3500HD

XL-4450



- High image quality
- Three preset FICE patterns available
- Applicable to full HD wide monitors

EPX-3500HD

Electronic Video Endoscopy System



High definition images and superior functions for advanced endoscopic diagnosis and treatment.

High image quality

VP-3500HD supports diagnoses with its advanced image processing technologies not incorporated in existing standard-class models. It offers clear images by using superior functions such as structure emphasis, color emphasis, automatic light control and anti-blur.


Specifications

Three preset FICE* patterns available



Three patterns of FICE*, which enhances the color tone of the endoscopic images by estimating spectral images of specific wavelength components from the original image by image processing, are preset and can be easily operated using the scope switch.

*Flexible spectral Imaging Color Enhancement (Image-Enhanced Endoscopy)

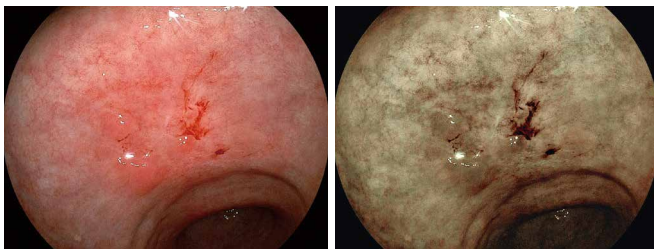

White light image

FICE 0

Applicable to full HD wide monitors

Captured image is displayed without overlapping with the real-time view when a full HD wide monitor is used.



Full HD monitor

5×4 monitor

XL-4450 Specification	
L	

FUJIFILM

FUJIFILM Corporation

26-30, NISHIAZABU 2-CHOME, MINATO-KU, TOKYO 106-8620, JAPAN
<http://www.fujifilm.com/>